

USAARL Report No. 92-6

DTIC
STANDARD
C

AD-A248 355



**Test and Evaluation Report
of the Physio Control Blood Pressure Monitor
Model LIFESTAT® 200**

By

Jeffrey D. Haun (Project Officer)

and

Joseph R. Licina (Project Officer)

Bill Olding (UES, Inc.)
Martin Quattlebaum (UES, Inc.)

Biodynamics Research Division

September 1991

92-09049



92 4 08 016

Approved for public release; distribution unlimited.

United States Army Aeromedical Research Laboratory
Fort Rucker, Alabama 36362-5292

Notice

Qualified requesters

Qualified requesters may obtain copies from the Defense Technical Information Center (DTIC), Cameron Station, Alexandria, Virginia 22314. Orders will be expedited if placed through the librarian or other person designated to request documents from DTIC.

Change of address

Organizations receiving reports from the U.S. Army Aeromedical Research Laboratory on automatic mailing lists should confirm correct address when corresponding about laboratory reports.


Disposition

Destroy this document when it is no longer needed. Do not return it to the originator.

Disclaimer


The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other official documentation. Citation of trade names in this report does not constitute an official Department of the Army endorsement or approval of the use of such commercial items.

Reviewed:

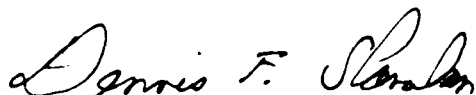


DENNIS F. SHANAHAN
LTC, MC, MFS
Director, Biodynamics
Research Division

Released for publication:

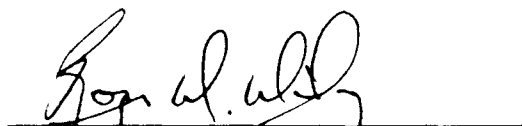

ROGER W. WILEY, O.D., Ph.D.
Chairman, Scientific
Review Committee



DAVID H. KARNEY
Colonel, MC, SFS
Commanding

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE

REPORT DOCUMENTATION PAGE

Form Approved
OMB No. 0704-0188

1a. REPORT SECURITY CLASSIFICATION UNCLASSIFIED			1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY			3. DISTRIBUTION/AVAILABILITY OF REPORT Approved for public release, distribution unlimited		
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE					
4. PERFORMING ORGANIZATION REPORT NUMBER(S) USAARL Report No. 92-6			5. MONITORING ORGANIZATION REPORT NUMBER(S)		
6a. NAME OF PERFORMING ORGANIZATION U.S. Army Aeromedical Research Laboratory		6b. OFFICE SYMBOL (If applicable) SGRD-UAD-IE		7a. NAME OF MONITORING ORGANIZATION U.S. Army Medical Research and Development Command	
6c. ADDRESS (City, State, and ZIP Code) P.O. Box 577 Fort Rucker, AL 36362-5292			7b. ADDRESS (City, State, and ZIP Code) Fort Detrick Frederick, MD 21702-5012		
8a. NAME OF FUNDING/SPONSORING ORGANIZATION		8b. OFFICE SYMBOL (If applicable)		9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER Partial effort under contract No. DAMD 17-86-C-6215	
9c. ADDRESS (City, State, and ZIP Code)			10. SOURCE OF FUNDING NUMBERS		
			PROGRAM ELEMENT NO. 0603807A	PROJECT NO. 3M463807D836	TASK NO. LC
11. TITLE (Include Security Classification) Test and Evaluation Report of the Physio Control Blood Pressure Monitor Model LIFESTAT® 200					
12. PERSONAL AUTHOR(S) Jeffrey D. H...n, Joseph R. Licina, Bill Olding, Martin Quattlebaum					
13a. TYPE OF REPORT Final		13b. TIME COVERED FROM _____ TO _____		14. DATE OF REPORT (Year, Month, Day) 1991 September	
15. PAGE COUNT 63					
16. SUPPLEMENTARY NOTATION					
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number) Electromagnetic compatibility, test and evaluation, aeromedical equipment		
FIELD	GROUP	SUB-GROUP			
19. ABSTRACT (Continue on reverse if necessary and identify by block number) The Physio Control Blood Pressure Monitor Model LIFESTAT® 200 was tested for electromagnetic interference/compatibility in the UH-60A helicopter under the U.S. Army Test and Evaluation of Aeromedical Equipment program. The tests were conducted using current military and industrial standards and procedures for electromagnetic interference/compatibility and human factors. The Physio Control Blood Pressure Monitor Model LIFESTAT® 200 was found to be compatible with U.S. Army medical evacuation UH-60A Blackhawk.					
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT <input type="checkbox"/> DTIC USERS				21. ABSTRACT SECURITY CLASSIFICATION UNCLASSIFIED	
22a. NAME OF RESPONSIBLE INDIVIDUAL Chief, Scientific Information Center				22b. TELEPHONE (Include Area Code) (205) 255-6907	
				22c. OFFICE SYMBOL SGRD-UAX-SI	

Table of contents

SECTION	PAGE
1. <u>EXECUTIVE DIGEST</u>	
1.1 Test objectives	1-1
1.2 Testing authority	1-2
1.3 Scope	1-2
1.4 Material description	1-3
1.5 Summary	1-3
1.6 Conclusion	1-5
2. <u>SUBTESTS</u>	
2.1 Initial inspection	2-1
2.2 Battery life evaluation	2-1
2.3 Electrical safety evaluation	2-2
2.4 Human factors evaluation (laboratory)	2-3
2.5 Altitude (low pressure) test	2-3
2.6 Vibration test	2-4
2.7 High temperature test	2-6
2.8 Low temperature test	2-7
2.9 Humidity test	2-8
2.10 Electromagnetic characteristics test	2-9
2.11 In-flight human factors evaluation	2-12
2.12 In-flight EMI/EMC characteristics test	2-13
3. <u>SUPPORTING DOCUMENTATION</u>	
3.1 Detailed test information	3-1

Table of contents (Continued)

3.2	Test data	3-3
3.3	Criteria, significant problems, and suggested improvements	3-36
3.4	References	3-39
3.5	Abbreviations	3-41
3.6	List of manufacturers.	3-43
3.7	Distribution list	3-44

Section 1. Executive digest

The Army program for Test and Evaluation of Aeromedical Equipment uses existing military standards (MIL-STD) and collective professional expertise to test and evaluate selected medical equipment proposed for use aboard Army aircraft. Equipment meeting these standards ensures the safety of the crew, patients, and aircraft by eliminating risks due to: (1) Interference by the medical equipment with aircraft systems/subsystems operation, (2) the aircraft system's interference with the operation of the medical equipment, (3) the medical equipment's susceptibility to environmental exposure, or (4) physical and/or functional incompatibility while in use on board selected rotary-wing aircraft. This program tests both developmental and nondevelopmental (off the shelf) medical equipment destined for use aboard Army medical evacuation aircraft.

1.1 TEST OBJECTIVES

1.1.1 To determine if the medical equipment is complete and operational per the manufacturer's operating instructions.

1.1.2 To ensure the electrical safety of the medical equipment.

1.1.3 To ensure the equipment will function as designed throughout the rated battery operation time.

1.1.4 To ensure the safety of the operator, the patient, and the aircrew.

1.1.5 To assess design considerations which potentially could contribute to an operator error.

1.1.6 To determine if the medical equipment can function as designed in a low pressure environment.

1.1.7 To determine the ability of the medical equipment to withstand the vibrational stresses expected in a rotary-wing flight environment without degradation or malfunction.

1.1.8 To determine the ability of the medical equipment to be stored and operated in a high temperature environment.

1.1.9 To determine the ability of the medical equipment to be stored and operated in a low temperature environment.

1.1.10 To determine the ability of the medical equipment to operate satisfactorily for short periods during exposure to high humidity conditions.

1.1.11 To assess the levels of electromagnetic emissions produced by the medical equipment within selected frequency ranges.

1.1.12 To assess the minimum electromagnetic susceptibility levels of the medical equipment within selected frequency ranges.

1.1.13 To assess the physical and/or functional compatibility of the medical equipment while in use on board the aircraft.

1.1.14 To assess the electromagnetic interference (EMI) and electromagnetic compatibility (EMC) characteristics of the medical equipment with the host aircraft and its installed systems.

1.2 TESTING AUTHORITY

Research and Technology Work Unit Summary, dated 5 October 1989. Project number 3M463807D836, titled, Army Program for Testing and Evaluation of Equipment for Aeromedical Operations.

1.3 SCOPE

1.3.1 This test was conducted at the United States Army Aeromedical Research Laboratory (USAARL), Cairns Army Airfield (CAAF), and designated test flight areas in and around Fort Rucker, Alabama.

1.3.2 The USAARL UH-60A aircraft, serial number 88-26069, with subsystems delineated in paragraph 3.2.2, was configured with the Physio Control blood pressure monitor, model LIFESTAT® 200 and used as the test aircraft for the in-flight evaluation. The in-flight evaluation required 2.5 flight hours.

1.3.3 Laboratory testing was accomplished at USAARL using government furnished equipment (GFE) by Universal Energy Systems, Inc. (UES), under contract No. DAMD 17-86-C-6215.

1.3.4 Prior to flight testing, the following tests were accomplished: Acceptance inspection, equipment training, electromagnetic compatibility, human factors and safety, environmental compatibility, and in-flight compatibility.

1.3.5 An airworthiness release (AWR) dated 12 July 1990 was received from the U.S. Army Aviation Systems Command (AVSCOM) prior to the in-flight testing of the LIFESTAT® 200.

1.4 MATERIAL DESCRIPTION

The Physio Control LIFESTAT® 200* is a portable blood pressure (BP) and pulse measurement device. It is used for noninvasive determination of systolic, diastolic and mean arterial pressures, and pulse rate. The unit is capable of taking manually initiated or automatic measurements at timed intervals of 1, 2, 3, 5, 10, 15, or 30 minutes. The operation is controlled by a microprocessor-based system. The light emitting diode (LED) displays present the systolic, diastolic, and mean pressures in addition to the pulse rate, cycle time, elapsed time, and alarm settings. Membrane switches on the front panel select manual or automatic mode, cuff inflation or deflation, alarm settings, and alarm silence. The paper strip chart recorder is controlled by two membrane switches. The strip chart records systolic, diastolic and mean pressures and pulse rate, patient identification, date, and time for each measurement. The instrument is powered by either internal battery or AC line power which is selected by a rocker switch on the front panel. The battery is a rechargeable, sealed lead-acid type. The front panel of the instrument contains the LED displays, membrane switch controls, chart recorder output, and the power switch. A Luer fitting is provided on the side of the unit to connect the blood pressure cuff. A volume control, AC power input, and "SYSTEMS OUTPUT" connector are provided on the rear panel.

1.5 SUMMARY

1.5.1 Laboratory testing

1.5.1.1 Battery Life Evaluation: The LIFESTAT® 200 was set to take automatic measurements at 2-minute intervals with the recorder on. In this mode, the BP monitor operated for an average of 39 minutes from a fully charged battery. This value agrees with the operator manual specification of 30 minutes operating time in this mode.

1.5.1.2 Electrical Safety Evaluation: All measurements were within acceptable limits. No unsafe qualities were found in the LIFESTAT® 200. The limits for currents and resistances were in accordance with (IAW) the National Association of Fire Prevention (NAFP) standards.

* See manufacturer's list

Accession For	
NTIS GRANT	<input checked="checked" type="checkbox"/>
DVCS TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution/	
Availability Codes	
Avail and/or	
Dist	Special
A-1	

1.5.1.3 Human Factors Evaluation: The LIFESTAT® 200 was found to be satisfactory in all major categories of the evaluation criteria with one exception. There are no externally accessible fuses, circuit breakers, or calibration points. The "T" adaptor, which is used during calibration checks, must be supplied by the user.

1.5.1.4 Environmental Tests: The LIFESTAT® 200 can be expected to perform in a variety of environmental conditions. Its performance was found to be satisfactory in all stages of the environmental testing. The requirements for environmental tests are established in MIL-STD-810D, methods 500.2 (altitude), 514.3 (vibration), 501.2 (high temperature), 502.2 (low temperature), and 507.2 (humidity).

1.5.1.5 Radiated Emissions Tests (RE02): The LIFESTAT® 200 may be unsatisfactory for use in certain EMI sensitive environments. Narrowband (NB) and broadband (BB) radiated emissions were detected in the test frequency ranges. Some narrowband and broadband emissions exceeded the test limits. Emission limits are set forth in MIL-STD-461A, Notice 4.

1.5.1.6 Radiated Susceptibility Test (RS03): The LIFESTAT® 200 was not found to be susceptible to radio frequency interference in the testing range and magnitude.

1.5.1.7 Conducted Emissions Test (CE01, CE02, and CE04): Narrowband signals were detected in the frequency range 13 to 53 MHz, with magnitudes 0.2 to 13.3 dB over specification limits. Broadband emissions were detected in the frequency range 15.6 to 53 MHz, with magnitudes 3 to 4.5 dB over specification limits.

1.5.1.8 Conducted Susceptibility Test (CS02 and CS06): Noise generated on the power lines by the LifeStat® 200 was greater than the test signal level. The LifeStat® 200 was not affected by the presence of the test spikes on its power lines.

1.5.2 In-flight testing

1.5.2.1 During the in-flight human factors evaluation, the LIFESTAT® 200 was found to be satisfactory in all categories of the evaluation criteria with one exception. There are no externally accessible fuses, circuit breakers, or calibration points.

1.5.2.2 The aircraft and its subsystems were not adversely affected by the operation of the LIFESTAT® 200 in any of the prescribed flight test modes.

1.5.2.3 The LIFESTAT® 200 was not affected by the aircraft and its subsystems during the in-flight testing.

1.6 CONCLUSIONS

Based on the results of laboratory and in-flight testing, the LIFESTAT® 200 was found to be compatible with the U.S. Army medical evacuation UH-60A Blackhawk with the subsystems listed in paragraph 3.2.2.

Section 2. Subtests

2.1 INITIAL INSPECTION

2.1.1 Objective

To determine if the LIFESTAT® 200 is complete and operational for testing per the manufacturer's operating instructions.

2.1.2 Criteria

2.1.2.1 The physical inventory is conducted solely for investigation and documentation.

2.1.2.2 The LIFESTAT® 200 will display consistent and accurate measurements as an acceptable performance test.

2.1.3 Test procedure

2.1.3.1 A complete physical inventory of the LIFESTAT® 200 was completed per the manufacturer's equipment list.

2.1.3.2 An operational validation test of the LIFESTAT® 200 was conducted per the manufacturer's operating instructions by USAARL's medical maintenance personnel.

2.1.4 Test findings

2.1.4.1 The LIFESTAT® 200 was inventoried and found to be complete.

2.1.4.2 The LIFESTAT® 200 operated as prescribed in the manufacturer's operating manual P/N 802611-02. Criteria met.

2.2 BATTERY LIFE EVALUATION (Laboratory)

2.2.1 Objective

To ensure the equipment will function as designed throughout the rated battery operation time.

2.2.2 Criterion

Verify manufacturer's specified full power internal battery life expectancy of 30 minutes during continuous operation in the 2-minute cycle mode, with the recorder on, and blood pressure measurements taken automatically at 2-minute intervals.

2.2.3 Test procedure

2.2.3.1 Charging and operation cycles were conducted in ambient room conditions.

2.2.3.2 The LIFESTAT® 200 was operated continuously using its fully charged internal battery in the 2-minute cycle mode until a low battery indication occurred. The depletion time was noted and the battery was recharged. This procedure was repeated three times.

2.2.4 Test findings

The test was conducted using the fully charged internal battery. The average operating time in testing was 39 minutes at room temperature. This exceeds manufacturer specification of 30 minutes. Criterion met.

2.3 ELECTRICAL SAFETY EVALUATION

2.3.1 Objective

To ensure the electrical safety, by evaluation of case-to-ground resistance and case-to-ground current leakage, of the LIFESTAT® 200.

2.3.2 Criterion

The LIFESTAT® 200 shall meet the standards established in National Association of Fire Prevention (NAFP) 99 for electrical safety of medical equipment.

2.3.3 Test procedure

Measurements in the electrical safety evaluation were made, with a Neurodyne-Dempsey model 431F electrical safety analyzer* IAW the procedures described in Technical Bulletin (TB) Number 38-750-2. Case-to-ground resistance and various case-to-ground leakage currents were measured. Leakage currents were measured using a 10 by 20 centimeter aluminum foil sheet taped flush to the equipment case. Checks were made for safety concerns such as case integrity, breaks in power cord insulation, and connectors.

2.3.4 Test findings

Grounding conductor resistance was 53.3 milliohms and maximum case leakage current was 12.2 microamperes. These measurements are below the limits specified in NAFP 99. Criterion met.

2.4 HUMAN FACTORS EVALUATION (Laboratory)

2.4.1 Objectives

2.4.1.1 To assure the safety of the operator, the potential patient, and the aircrew.

2.4.1.2 To assess the design considerations which could potentially contribute to an operator error.

2.4.2 Criterion

The LIFESTAT® 200 must be rated satisfactory in all major categories of the evaluation. These include: (1) visual displays, (2) controls, (3) maintainability, (4) conductors, (5) fasteners, (6) test points, (7) test equipment, (8) fuses and circuit breakers, (9) labels and coding, and (10) safety.

2.4.3 Test procedure

2.4.3.1 The evaluation was conducted in a laboratory under fluorescent lighting and ambient room conditions.

2.4.3.2 The LIFESTAT® 200 was operated according to prescribed instructions through its full range of functions.

2.4.4 Test finding

The LIFESTAT® 200 was found to be satisfactory in all of the evaluation criteria with one exception. There are no externally accessible fuses, circuit breakers, or calibration points. Criterion partially met.

2.5 ALTITUDE (LOW PRESSURE) TEST [IAW MIL-STD-810D, METHOD 500.2]

2.5.1 Objective

To determine if the LIFESTAT® 200 can function as designed in a low pressure environment.

2.5.2 Criterion

The LIFESTAT® 200 will display consistent and accurate measurements while exposed to an altitude equivalency of 15,000 feet above sea level.

2.5.3 Test procedure

2.5.3.1 A pretest performance check was conducted to ensure proper operation of the LIFESTAT® 200.

2.5.3.2 The altitude test was performed in a Tenney Engineering model 64S altitude chamber. This test is based on MIL-STD-810D, Method 500.2. The LIFESTAT® 200 was turned on in the standby mode and placed on the floor of the chamber. Chamber pressure was decreased to 420 mmHg (15,000 ft equivalent altitude) over a 15-minute period, held constant for 60 minutes, then raised, at 1500 fpm, to ambient conditions (760 mm Hg) over a 10-minute period. There are no provisions for the control of temperature or humidity inside this chamber.

2.5.3.3 A posttest performance check was conducted to ensure proper operation of the LIFESTAT® 200 after the exposure to low pressure.

2.5.4 Test findings

2.5.4.1 The pretest performance check met criterion 2.1.2.2.

2.5.4.2 No failures in the performance of the LIFESTAT® 200 were noted before, during, or after the altitude test. Criterion met.

2.5.4.3 The posttest performance check met criterion 2.1.2.2.

2.6 VIBRATION TEST [IAW MIL-STD-810D, METHOD 514.3]

2.6.1 Objective

To determine the ability of the LIFESTAT® 200 to withstand the vibrational stresses expected in a rotary-wing environment without degradation or malfunction.

2.6.2 Criterion

The LIFESTAT® 200 will remain operational and be able to display consistent and accurate measurements while exposed to vibrational stresses.

2.6.3 Test procedure

2.6.3.1 A pretest performance check was conducted to ensure proper operation of the LIFESTAT® 200.

2.6.3.2 The vibration test was performed using an Unholtz-Dickey model TA115-40/CSTA vibration test system*. It is a single-axis system with an electromagnetic driver unit. The test consisted of sinusoidal vibrations superimposed on random vibrations over a frequency range of 500 Hz, as shown below. These vibrations are derived from measurements taken on the floor under the copilot's seat in a UH-1 helicopter traveling at 120 knots. The reference spectrum breakpoints are from MIL-STD-810D, Method 514.3; reference spectrum levels are based on field measurements with a conservatism factor of 1.5. Independent tests were conducted in the X, Y, and Z axes.

Z-axis

duration: 60 minutes
broadband intensity: $0.4506 G_{rms}$
random vibration: initial slope : 99.00 dB/Hz
5 Hz level: $0.00006210 G_{sqr/Hz}$
100 Hz level: $0.0006210 G_{sqr/Hz}$
300 Hz level: $0.0006210 G_{sqr/Hz}$
500 Hz level: $0.00006210 G_{sqr/Hz}$
final slope: -99.00 dB/oct
sinusoidal vibration: $.5450 G_{pk}$ at 11.25 Hz
 $.1690 G_{pk}$ at 22.50 Hz
 $.1200 G_{pk}$ at 33.75 Hz
 $.0310 G_{pk}$ at 45.00 Hz
 $.0530 G_{pk}$ at 56.25 Hz

X and Y axes

duration: 60 minutes each
broadband intensity: $0.3099 G_{rms}$
random vibration: initial slope: 99.00 dB/oct
5 Hz level: $0.00002920 G_{sqr/Hz}$
100 Hz level: $0.0002920 G_{sqr/Hz}$
300 Hz level: $0.0002920 G_{sqr/Hz}$
500 Hz level: $0.00002920 G_{sqr/Hz}$
final slope: -99.00 dB/oct
sinusoidal vibration: $.3200 G_{pk}$ at 11.25 Hz
 $.0670 G_{pk}$ at 22.50 Hz
 $.0950 G_{pk}$ at 33.75 Hz
 $.0350 G_{pk}$ at 45.00 Hz
 $.0770 G_{pk}$ at 56.25 Hz

The LIFESTAT® 200 was strapped to the vibration table fixture, and its performance was evaluated before, during, and after exposure to vibration.

2.6.3.3 A posttest performance check was conducted to ensure proper operation of the LIFESTAT® 200.

2.6.4 Test findings

2.6.4.1 The pretest performance check met criterion 2.1.2.2.

2.6.4.2 No failures in the performance of the LIFESTAT® 200 occurred before, during, or after exposure to vibration. Criterion met.

2.6.4.3 The posttest performance check met criterion 2.1.2.2.

2.7 HIGH TEMPERATURE TEST [IAW MIL-STD-810D, METHOD 501.2]

2.7.1 Objective

To determine the ability of the LIFESTAT® 200 to be stored and operated in a high temperature environment.

2.7.2 Criteria

2.7.2.1 The LIFESTAT® 200 will display consistent and accurate measurements during the high temperature operation check.

2.7.2.2 The LIFESTAT® 200 will display consistent and accurate measurements after the high temperature storage cycle.

2.7.3 Test procedure

2.7.3.1 A pretest performance check was conducted to ensure proper operation of the LIFESTAT® 200.

2.7.3.2 The high temperature test was conducted in a Tenney Engineering model ZWUL-10107D walk-in controlled environment chamber. This test is based on MIL-STD-810D, Method 501.2. For the high temperature operation test, the LIFESTAT® 200 was turned on in the standby mode and placed on the floor of the environmental chamber. The chamber temperature was raised to 49°C and the humidity was stabilized at a maximum of 20 percent RH within 15 minutes. The environmental control system is capable of regulating temperature within $\pm 2^\circ\text{C}$ and humidity within ± 5 percent RH. Temperature and humidity were held constant for 2 hours. At 30-minute intervals, the chamber door was opened briefly to minimize the change in chamber conditions during performance checks. After the operational test, the LIFESTAT® 200 was allowed to return to ambient conditions over a 30-minute period.

2.7.3.3 A posttest performance check was conducted to ensure proper operation of the LIFESTAT® 200.

2.7.3.4 The LIFESTAT® 200 was stored (not operated) at temperatures of 63°C for 1 hour, 71°C for 4 hours, then again at 63°C for 1 hour. The chamber and LIFESTAT® 200 then were returned to ambient conditions over a 30-minute period.

2.7.3.5 A post storage performance check was conducted to ensure proper performance of the LIFESTAT® 200.

2.7.4 Test findings

2.7.4.1 The pretest performance check met criterion 2.1.2.2.

2.7.4.2 No operational failures occurred during the high temperature test. Criterion met.

2.7.4.3 The posttest performance check met criterion 2.1.2.2.

2.7.4.4 The LIFESTAT® 200 functioned properly after the high temperature storage test. Criterion met.

2.8 LOW TEMPERATURE TEST [IAW MIL-STD-810D, METHOD 502.2]

2.8.1 Objective

To determine the ability of the LIFESTAT® 200 to be stored and operated in a low temperature environment.

2.8.2 Criteria

2.8.2.1 The LIFESTAT® 200 will display consistent and accurate measurements during the low temperature operation check.

2.8.2.2 The LIFESTAT® 200 will display consistent and accurate measurements after the low temperature storage cycle.

2.8.3 Test procedure

2.8.3.1 A pretest performance check was conducted to ensure proper operation of the LIFESTAT® 200.

2.8.3.2 The LIFESTAT® 200 was placed on the floor of the environmental chamber and the temperature was lowered to 0°C within 25 minutes. The environmental control system is capable of regulating temperature within 2°C. Humidity cannot be controlled in the chamber at freezing temperatures. The temperature was held constant for 2 hours. The chamber door was opened briefly every 30 minutes to minimize the change in chamber conditions and

a performance check was conducted. The chamber temperature then was raised to ambient temperature within a 30-minute period.

2.8.3.3 A posttest performance check was conducted to ensure proper operation of the LIFESTAT® 200.

2.8.3.4 The LIFESTAT® 200 was "stored" in a nonoperational mode. The LIFESTAT® 200 was placed on the floor of the environmental test chamber and the temperature was lowered to -46°C for 6 hours. The chamber then was raised to ambient temperature over a 30-minute period.

2.8.3.5 A post storage performance check was conducted to ensure proper operation of the LIFESTAT® 200.

2.8.4 Test findings

2.8.4.1 The pretest performance check met criterion 2.1.2.2.

2.8.4.2 No operational failures occurred during the low temperature test. Criterion met.

2.8.4.3 The posttest performance check met criterion 2.1.2.2.

2.8.4.4 The LIFESTAT® 200 functioned properly after the low temperature storage test. Criterion met.

2.9 HUMIDITY TEST [IAW MIL-STD-810D, METHOD 507.2]

2.9.1 Objective

To determine the ability of the LIFESTAT® 200 to operate satisfactorily for short periods of time during exposure to highly humid conditions.

2.9.2 Criterion

The LIFESTAT® 200 will display consistent and accurate measurements while exposed to a high humidity environment.

2.9.3 Test procedure

2.9.3.1 A pretest performance check was conducted to ensure the proper operation of the LIFESTAT® 200.

2.9.3.2 The humidity test was conducted in a Tenney Engineering model ZWUL-10107D walk-in controlled environment chamber*. This test is based on MIL-STD-810D, Method 507.2. For the humidity test, the LIFESTAT® 200 was placed ready for operation on the floor of the environmental chamber. The chamber temperature was

raised to a temperature of 30°C and a relative humidity of 95 percent within 25 minutes. Temperature and relative humidity were maintained for 4 hours. The environmental control system is capable of regulating temperature within $\pm 2^\circ\text{C}$ and humidity within ± 5 percent RH. At 45-minute intervals the performance of the blood pressure monitor was checked. The chamber door was opened briefly to minimize the change in chamber conditions. The chamber and the LIFESTAT® 200 were returned to ambient conditions before the posttest performance validation check was conducted.

2.9.3.3 A posttest performance check was conducted to ensure the proper operation of the LIFESTAT® 200.

2.9.4 Test findings

2.9.4.1 The pretest performance check met criterion 2.1.2.2.

2.9.4.2 No failures were noted in the LIFESTAT® 200 performance checks conducted during the exposure to the high humidity environment. Criterion met.

2.9.4.3 The posttest performance check met criterion 2.1.2.2.

2.10 ELECTROMAGNETIC CHARACTERISTICS TEST [IAW MIL-STD-461A, Notice 4, MIL-STD-462, Notice 3, and MIL-STD-704C]

2.10.1 Objectives

2.10.1.1 To assess the maximum levels of radiated electromagnetic emissions produced by the LIFESTAT® 200 in the 14 kHz to 12.4 GHz frequency range.

2.10.1.2 To assess the tolerances of radiated electromagnetic susceptibility of the LIFESTAT® 200 within the 10 kHz to 10 GHz electric field.

2.10.1.3 To assess the maximum levels of conducted electromagnetic emissions produced by the LIFESTAT® 200 in the 10 kHz to 50 MHz frequency ranges.

2.10.1.4 To assess the tolerances of conducted electromagnetic susceptibility of the LIFESTAT® 200 within the range of 50 kHz to 400 MHz and power spikes.

2.10.2 Criteria

2.10.2.1 The LIFESTAT® 200 will not produce emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraph 6.13.

2.10.2.2 The LIFESTAT® 200 will not malfunction when it is subjected to radiated emissions as specified in MIL-STD-461A, Notice 4, paragraph 6.20.

2.10.2.3 The LIFESTAT® 200 shall not conduct emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraphs 6.1 and 6.2.

2.10.2.4 The LIFESTAT® 200 shall not malfunction when it is subjected to conducted emissions as specified in MIL-STD-461A, Notice 4, paragraphs 6.7 and 6.10.

2.10.3 Test procedure

2.10.3.1 The radiated emissions test was performed according to MIL-STD-462, Notice 3, Method RE02. The LIFESTAT® 200 was positioned on a wooden test stand inside the EMI chamber, 1 meter away from the receiving antennas. The antennas were mounted for both vertical and horizontal polarities and connected to EMI receivers. The LIFESTAT® 200 was connected through an extended tube to a cuff outside the chamber. The cuff was placed around a test engineer's arm while the LIFESTAT® 200 took blood pressure measurements at 2-minute intervals. While the LIFESTAT® 200 was operating, the frequency spectrum (14 kHz to 12.4 GHz) was scanned for emissions. The LIFESTAT® 200 was operated with both AC and battery power.

2.10.3.2 The radiated susceptibility test was performed according to MIL-STD-462, Notice 3, Method RS03. The LIFESTAT® 200 was positioned on a wooden test stand inside the EMI chamber 1 meter away from the transmitting antennas. The antennas were mounted for both vertical and horizontal polarities and connected to radio frequency (RF) transmitters. The LIFESTAT® 200 was connected through an extended tube to a cuff outside the chamber. The cuff was placed around a test engineer's arm while the LIFESTAT® 200 took blood pressure measurements at 2-minute intervals. While the LIFESTAT® 200 was operating, it was monitored for faulty operation during exposures to fields of 1 V/m from 10 kHz to 2 MHz, and 5 V/m from 2 to 30 MHz, 10 V/m from 30 MHz to 2 GHz, and 5V/m from 2 to 10 GHz. The LIFESTAT® 200 was operated with AC power only.

2.10.3.3 The conducted emissions tests were performed according to MIL-STD-462, Notice 3, Methods CE02 and CE04. The LIFESTAT® 200 was placed on a grounded, copper covered workbench. The top of the workbench was 1 meter from floor level, 1.37 meters long and 0.81 meters wide. Power was supplied via a pair of line impedance stabilization networks (LISN's) and a test jig. The test jig is a wooden tray with two power receptacles and two slots to hold current probes in place around power supply conductors. While the LIFESTAT® 200 was operating, the frequency range

(10 kHz to 50 MHz) was scanned for emissions conducted in the power cable from the LIFESTAT® 200.

2.10.3.4 The conducted susceptibility spike test was performed according to MIL-STD-462, Notice 3, Method CS06, on a chemical resistant counter top. Power was supplied via a customized metal connection box. The connection box has two power receptacles and four banana jacks on its front panel. Connections to the individual power lines are made in series through the banana jacks. Transient spikes of 100 volts, 10 microseconds were generated with a Solar Electronics model 8282-1 transient pulse generator* and induced onto the power leads at the connection box banana jacks. The spikes were monitored with a Tektronix 2235 oscilloscope* connected to a power receptacle on the connection box. The LIFESTAT® 200 was plugged into the other receptacle on the connection box and placed in operation. It was observed for correct operation and visual displays while it was subjected to the power line spikes.

2.10.3.5 The conducted susceptibility test was performed according to MIL-STD-462, Notice 3, Method CS02. The LIFESTAT® 200 was placed on a grounded, copper covered workbench. Radio frequency interference was induced on the power leads and measured at the LIFESTAT® 200 power cable. The frequency of the interference was incremented over the 50 kHz to 400 MHz range while the LIFESTAT® 200 was operated. It was observed for correct operation and visual displays while it was subjected to the radio interference on the power leads. Each frequency was held for 15 seconds.

2.10.4 Test findings

2.10.4.1 During the radiated emissions test, emissions which exceeded specification limits of MIL-STD-461, Notice 4, were detected in the narrowband frequency range 2 to 255 MHz, with magnitudes 0.3 dB to 41 dB over the specification limits, and in the broadband frequency range 15 to 45 MHz, with magnitudes 1.8 dB to 21.3 dB over the specification limits. Criterion partially met.

2.10.4.2 The LIFESTAT® 200 was not found to be susceptible to radio frequency interference in the testing range and magnitude. Criterion met.

2.10.4.3 Narrowband signals were detected in the frequency range 13 to 53 MHz, with magnitudes 0.2 to 13.3 dB over specification limits. Broadband emissions were detected in the frequency range 15.6 to 53 MHz, with magnitudes 3 to 4.5 dB over specification limits. Criterion partially met.

2.10.4.4 Noise generated on the powerlines by the LifeStat® 200 was greater than the test signal level. The LifeStat® 200 was not susceptible to test spikes during the conducted susceptibility tests. Criterion partially met.

2.11 IN-FLIGHT HUMAN FACTORS EVALUATION

2.11.1 Objective

To assess the physical and/or functional compatibility of the LIFESTAT® 200 while in use onboard the aircraft.

2.11.2 Criterion

The flight surgeon will be able to operate the LIFESTAT® 200 without physical or functional restrictions aboard the aircraft. Major areas of concern include: Proper operation, visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.

2.11.3 Test procedure

2.11.3.1 A human factors evaluation was performed IAW MIL-STD-1472D, AAMI Human Factors Engineering Guidelines, and UL-544 to ensure the compatibility of the LIFESTAT® 200 and the in-flight environment. The flight surgeon conducted the test wearing a flight suit, flight gloves, and an SPH-4 flight helmet. An evaluation of the compatibility with the nuclear, biological, and chemical (NBC) protective equipment was not conducted. Due to restrictions of the AWR, testing was conducted during daylight hours only.

2.11.3.2 The LIFESTAT® 200 was placed on the floor of the aircraft next to the bottom pan of the litter carousel which was configured for four patients. The litter carousel was flown in the "load" position (perpendicular to the long axis of the helicopter). The LIFESTAT® 200 was tested with the cuff placed on the right arm of a simulated patient laying in the bottom pan of the litter carousel. The LIFESTAT® 200 was tested in both the defibrillation and monitoring modes in the flight scenarios noted in section 3.2).

2.11.4 Test findings

During the in-flight human factors evaluation, the LIFESTAT® 200 was found to be satisfactory in all but one of the categories of the evaluation criteria. A deficiency was noted in the lack of externally accessible fuses, circuit breakers, or calibration points as was noted in the laboratory evaluation (paragraph 1.5.1.3). The LED display was readable in all but direct sunlight which was not considered a problem in using the device. Another concern is the inability to vary the display intensity which could present a problem if used with night vision devices. Criterion partially met.

2.12 IN-FLIGHT EMI/EMC CHARACTERISTICS TEST

2.12.1 Objective

To assess the EMI/EMC characteristics of the LIFESTAT® 200 with the host aircraft and its installed systems.

2.12.2 Criteria

2.12.2.1 The LIFESTAT® 200 will not radiate EMI to disrupt or interfere with other equipment or systems aboard the aircraft.

2.12.2.2 The aircraft will not radiate EMI to disrupt or interfere with the LIFESTAT® 200's operation.

2.12.3 Test procedure

A qualitative EMI/EMC assessment was performed with both the LIFESTAT® 200 and the aircraft operating as source and victim. The LIFESTAT® 200 and applicable aircraft instruments and systems were monitored for unusual operation, readings, surges, or power anomalies for each checklist item (see pages 3-5 through 3-8).

2.12.4 Test findings

2.12.4.1 There were no adverse instances of EMI/EMC noted with the LIFESTAT® 200 acting as either the source or victim. Criterion met.

2.12.4.2 There were no adverse instances of EMI/EMC noted with the aircraft acting as either the source or victim. Criterion met.

SECTION 3. Supporting documentation

3.1 DETAILED TEST INFORMATION

3.1.1 General information

3.1.1.1 LIFESTAT® 200 testing is not considered a major action significantly affecting the quality of the human environment and therefore qualifies for categorical exclusion A-28, appendix A, AR 200-1.

3.1.1.2 A safety pilot will be designated for each flight. Flight operations will be conducted IAW the aircraft operator's manual, appropriate aircrew training manuals, and test item technical data.

3.1.2 Material description

3.1.2.1 The Physio Control LIFESTAT® 200 is a portable blood pressure and pulse measurement device. It is used for noninvasive determination of systolic, diastolic, and mean arterial pressures, and pulse rate. The unit is capable of taking manually initiated or automatic measurements at timed intervals of 1, 2, 3, 5, 10, 15, or 30 minutes. The operation is controlled by a microprocessor-based system. LED displays present the systolic, diastolic, and mean pressures in addition to the pulse rate, cycle time, elapsed time, and alarm settings. Membrane switches on the front panel select manual or automatic mode, cuff inflation or deflation, alarm settings, and alarm silence. The paper strip chart recorder is controlled by two membrane switches. The strip chart records systolic, diastolic, and mean pressures and pulse rate, patient identification, date, and time for each measurement. The instrument is powered by either internal battery or AC line power which is selected by a rocker switch on the front panel. The battery is a rechargeable, sealed lead-acid type. The front panel of the instrument contains the LED displays, membrane switch controls, chart recorder output, and the power switch. A Luer fitting is provided on the side of the unit to connect the blood pressure cuff. A volume control, AC power input, and "SYSTEMS OUTPUT" connector are provided on the rear panel.

3.1.2.2 Method of Operation: The LIFESTAT® 200's operation is based on the oscillometric technique. Arterial pulsations acting against the inflated cuff are used to determine blood pressure and pulse rate. These pulsations are analyzed by an internal computer which determines systolic, diastolic and mean arterial pressures, and pulse rate. It is able to distinguish between real heart beats and a certain amount of motion artifact. With too much artifact the computer displays "---" rather than

incorrect information. Cuff inflation pressure is user variable to 165 mmHg (LO), or 220 mmHg (HI), or 50 to 290 mmHg (OVERRIDE) in the manual mode. Cuff inflation pressure is user variable to 165 mmHg (LO), or 220 mmHg (HI) in the automatic mode, with timed measurements at 1, 2, 3, 5, 10, 15, and 30-minute intervals.

3.1.2.3 Dimensions: 8.3 x 5.8 x 11.5 in (21.1 x 14.7 x 29.2 cm).

3.1.2.4 Weight: 15 lb (6.8 kg), not including accessories.


3.1.2.5 Standard accessories: Printer paper, writing pens, pediatric cuff, adult cuff, large adult cuff, thigh cuff, latex extension sets, operating and service manual operating instructions.

3.1.2.6 Power requirements: 90 to 132 VAC, 50 or 60 Hz, 70 watts maximum. Internal battery is a sealed lead-acid type, 9.8 V nominal, 2.5 Ah. Battery capacity is 30 minutes with a 16-hour full charge, in 2-minute cycle using adult cuff, 6 ft tube with printer on. Power cord is type SJTW-A, 16/3C AWG conductor, E58188 VW-1.

3.1.2.7 Environmental considerations: atmospheric pressure, less than 11000 ft; operating temperature, 5 to 45 degrees C; storage temperature, -30 to 65 degrees C; relative humidity, 0 to 95 percent.

3.2 TEST DATA

3.2.1 Photographic Description

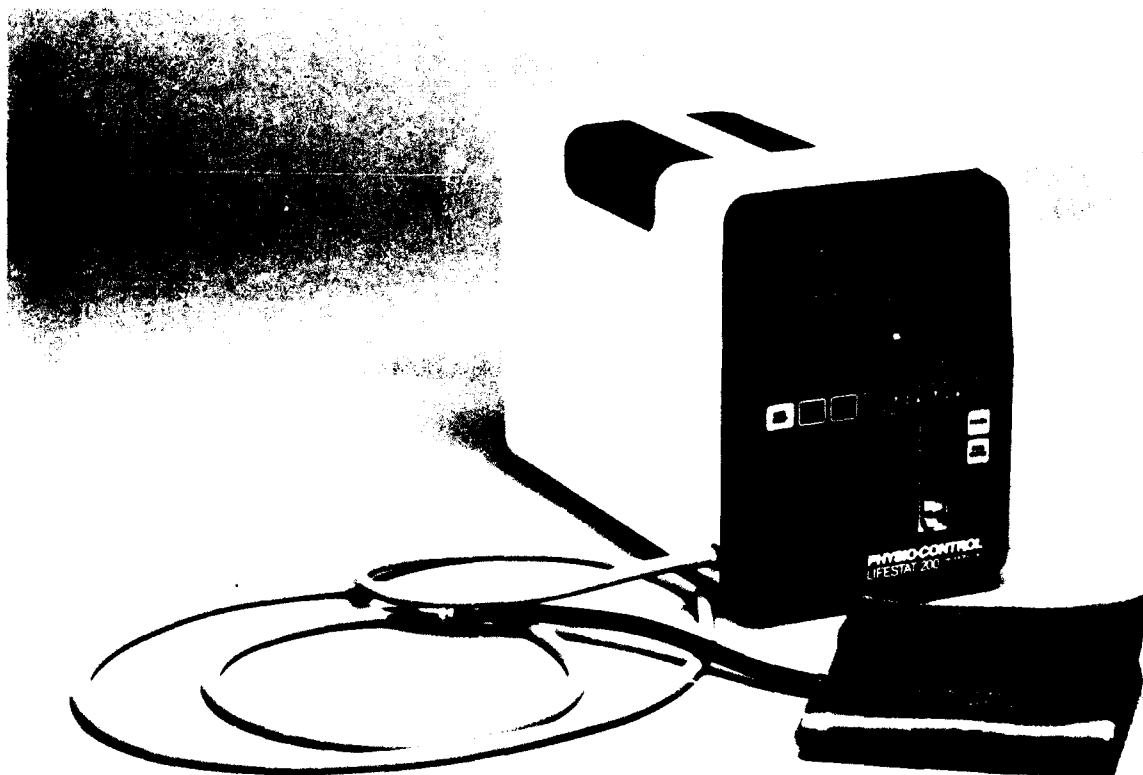

3.2.2 Aircraft equipment list

Item No.	Nomenclature
1	Receiver radio -- R-1496A/ARN-89 (automatic direction finder)
2	Displacement gyro -- CN-1314/A
3	Gyro directional -- CN-998/ASN-43
4	Signal data converter -- CV-3338/ASN-128
5	Receiver -- R-2139/ARN-143 (VOR/LOC/MB/GS)
6	Command instrument system processor -- 70600-01038-101
7	SAS amplifier -- 70901-02908-104 (flight control stability augmentation system)
8	Rate gyro -- TRU-2A/A
9	Amplifier, impedance -- AM-4859A/ARN-89
10	Cargo hook -- FE-7590-145
11	Receiver, radar -- RT-1193/ASN-128 (doppler navigation receiver)
12	Barometric altimeter -- AAU-31/A-1
13	Barometric altimeter -- AAU-32A
14	Receiver/transmitter -- RT-1300/ARC-186 (VHF-AM and/or FM radio)
15	UHF-AM radio set -- RT-1518/ARC-164
16	Interphone control -- C6533/ARC (aircraft intercom control)
17	Receiver/transmitter -- RT-1115D/APN-209 (radar altimeter)
18	Indicator altimeter -- ID-1917C/APN-209 (radar altimeter)
19	Control radio set -- C-7392A/ARN-89 (automatic direction finder)
20	Comparator signal data -- CM-482/ARC-186 (comparator for ARC-186)
21	Receiver/transmitter -- RT-1296A/APX-100 (transponder with IFF)
22	Computer display unit -- CP-1252/ASN-128 (doppler navigation system)
23	Compass set controller -- C-8021E/ASN75
24	Magnetic compass - standby -- MS-17983-4

3.2.3 In-flight test data card

DATA CARD FORMAT

GUIDELINE FOR DATA COLLECTION

IN-FLIGHT SUITABILITY TEST OF MEDICAL ITEMS

1. Installation/removal.	Suitable		Comments
	Yes	No	
a. Weight and balance (DD Form 365-4, Clearance Form F).	X		
b. Space/area allocation.			
(1) Operational requirements.	X		
(2) Storage requirements.	X		
c. Interface connections (safe, positive, secure).	X		
d. Installation/removal (expedient/easily achieved).	X		
e. Mounting/final config- uration (functional/stable).	X		
2. Operations and performance.	Suitable		Comments
	Yes	No	
a. Manufacturer's operating instruction.	X		
b. Medical item operation before aircraft run-up.	X		
c. System interface during aircraft engine run-up and medical item operation (EMI switchology checklist).	X		
(1) Aircraft voltage output.	X		

	Suitable		Comments
	Yes	No	
(2) Flight control function (UH-60).			
(3) Stabilator function (UH-60).	X		
(4) Radio communication vs medical item operation.			
(a) FM	X		
(b) UHF	X		
(c) VHF	X		
(5) Navigation equipment vs medical item operation.			
(a) Transponder	X		
(b) ADF	X		
(c) VOR	X		
(d) DOPPLER	X		
(6) Radar altimeter operation vs medical item operation.	X		

d. System interface during aircraft hover and medical item operation (EMI switchology checklist).

(1) Voltage output.	n/a
(2) Radio communication vs medical item operation.	
(a) FM	X
(b) UHF	X
(c) VHF	X

(3) Navigation equipment operation vs medical item operation.	Suitable		Comments
	Yes	No	
(a) Transponder	X		
(b) ADF	X		
(c) VOR	X		
(d) DOPPLER	X		
e. Flight mission profile vs medical item operation (EMI switchology checklist).			
(1) Straight and level (1000 ft MSL for 20 minutes).			
(a) Compatibility of flight mode and medical item operation.	X		
(b) Radio communication vs medical item opera- tion.			
<u>a.</u> FM	X		
<u>b.</u> UHF	X		
<u>c.</u> VHF	X		
(2) NOE (20 minutes). compatibility of flight mode and medical item operation.	X		
(3) FM homing (10 minutes).	X		
(4) DOPPLER navigation vs medical item operation.			
(a) Initialize function.	X		
(b) Fix function.	X		
(c) Update function.	X		

	Suitable Yes No	Comments
(5) VOR navigation vs medical item operation.	X	
(6) ILS approach vs medical item operation.	X	
f. Medical item operation after engine shutdown (external power source).	X	
g. Restrictions to the medical item's use (i.e., electrical connectors).	X	
h. Deviations from the labor- atory test results.		
(1) Electrical/ electronic.	None	
(2) Mechanical environment.	None	
(3) Human factors (user interface, controls, markings, lighting, egress).	None	
(4) Safety.	None	
3. Deviations from the in-flight test protocol.		
a. The VOR navigation portion of the in-flight test con- ducted at 2000 feet MSL due to air traffic control clearance.		

3.2.4 EMI switchology checklist

EMI SWITCHOLOGY CHECKLIST UH-60 AIRCRAFT

IN-FLIGHT SUITABILITY OF MEDICAL ITEMS

ENG INSTRUMENTS/CDU	No EMI Affect	EMI Affected Gnd Flt	Explanation
Fuel quantity	X		
Fuel indicator test	X		
XMSN oil temperature	X		
XMSN oil pressure	X		
#1 engine oil temperature	X		
#2 engine oil temperature	X		
#1 engine oil pressure	X		
#2 engine oil pressure	X		
#1 TGT	X		
#2 TGT	X		
#1 Ng speed	X		
#2 Ng speed	X		
CDU digits on/off	X		
CDU instruments dim	X		
ENG INSTRUMENTS/PLT PDU	No EMI Affect	EMI Affected Gnd Flt	Explanation
#1 engine RPM	X		
#2 engine RPM	X		
Rotor RPM	X		
#1 torque	X		
#2 torque	X		
ENG INSTRUMENTS/COPLT PDU	No EMI Affect	EMI Affected Gnd Flt	Explanation
#1 engine RPM	X		
#2 engine RPM	X		
Rotor RPM	X		
#1 torque	X		
#2 torque	X		

ENG CONTROLS	No EMI Affect	EMI Affected Gnd Flt	Explanation
#1 overspeed	X		
#2 overspeed	X		
RPM switch	X		
#1 engine anti-ice	X		
#2 engine anti-ice	X		
#1 inlet anti-ice	X		
#2 inlet anti-ice	X		

RADIO EQUIPMENT	No EMI Affect	EMI Affected Gnd Flt	Explanation
ICS, C-6533 ARC	X		
VHF-FM, ARC-186/115	X		
VHF-AM, ARC-186/115	X		
UHF-AM, ARC-164(V)	X		
Crypto, KY-28	Not installed		
Radio retransmissions PLN	Not installed		
Transponder, APX-100(V)	X		
KIT-1A/TSEC IFF computer	Not keyed with code		

MISSION EQUIPMENT	No EMI Affect	EMI Affected Gnd Flt	Explanation
RWR, APR-39(V)	Not installed		
IR CM, ALQ-144	Not installed		
Chaff dispenser, M-130	Not installed		
Cargo hook system	X		

HYDRAULIC CONTROL SYSTEM	No EMI Affect	EMI Affected Gnd Flt	Explanation
Backup hydraulic pump	X		
Servo off 1st stage/PLT	X		
Servo off 2nd stage/PLT	X		
Servo off 1st stage/COPLT	X		
Servo off 2nd stage/COPLT	X		
Hydraulic leak test	X		
Tail servo	X		
Boost servos	X		

FUEL SYSTEM	No EMI Affect	EMI Affected Gnd Flt	Explanation
Fuel pump switch	X		
Fuel boost pump #1	X		
Fuel boost pump #2	X		
Fuel cont panel ESSS	Not installed		

WARNING SYSTEM	No EMI Affect	EMI Affected Gnd Flt	Explanation
Low rotor RPM	X		
Master caution	X		
Caution advisory	X		
Fire warning	X		
AFCS	X		
Stabilator	X		
#1 engine out	X		
#2 engine out	X		

NAVIGATION INSTRUMENTS	No EMI Affect	EMI Affected Gnd Flt	Explanation
ADF	X		
Magnetic compass	X		
CONUS NAV, ARN-123	X		
DOPPLER, ASN-128	X		
Gyro mag compass (PLT)	X		
Gyro mag compass (COPLT)	X		
Compass cont panel, ASN-75	X		
HSI	X		

FLIGHT INSTRUMENTS	No EMI Affect	EMI Affected Gnd Flt	Explanation
Radar altimeter	X		
Stabilator pos indicator	X		
VSI	X		
CIS mode select	X		
SAS 1	X		
SAS 2	X		
FPS	X		
Trim	X		
Go-around enable	X		
Cyclic trim release	X		
Cyclic stick trim	X		
ALR encoder	X		

FLIGHT INSTRUMENTS (CONT)	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
HSI/VSI mode select (PLT)				
DPLR	X			
VOR/ILS	X			
BACK CRS	X			
FM HOME	X			
TURN RATE	X			
CRS HDG	X			
VERT GYRO	X			
BRG 2	X			
HSI/VSI Mode Select (COPLT)				
DPLR	X			
VOR/ILS	X			
BACK CRS	X			
FM HOME	X			
TURN RATE	X			
CRS HDG	X			
VERT GYRO	X			
BRG 2	X			
MISCELLANEOUS EQUIPMENT	No EMI Affect	EMI Affected Gnd Flt		Explanation
Blade deice	Not tested			Ambient tempera- ture was out of test lim- its.
Windshield anti-ice	X			
Pitot heat	X			
Vent blower	X			
Windshield wiper	X			
Heater	X			
APU	X			
Generator #1	X			
Generator #2	X			
Generator APU	X			
Air source heat start	X			
Tail wheel lock	X			
Gyro erect	X			

LIGHTING	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
Cockpit utility	X			
Cockpit flood	X			
Cabin dome	X			
Search light	X			
Search light control	X			
Landing light	X			
Flt instr lights (PLT)	X			
Flt instr lights (COPLT)	X			
Nonflight instr lights	X			
Console lights, upper	X			
Console lights, lower	X			
Position lights	X			
Formation lights	X			
Anticollision lights	X			
NVG lighting	X			

3.2.5 Battery life evaluation

Battery Life Evaluation Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 00003850
Military item number: None

Options installed: None

Manufacturer battery life specification: Approximately 30 minutes in 2-minute cycle, recorder on, when fully charged.

Specified battery recharge time: 16 hours to fully charge depleted battery; 2 hours to provide 70 percent of full charge capacity.

Specified mode of operation under battery power: 2-minute cycle mode, in which automatic blood pressure measurements are taken at 2-minute intervals and recorder on.

Overall performance: Pass

Measurements: The unit averaged 39 minutes of operation.

Comments: The unit was operated continuously in the 2-minute cycle mode until a low battery indication occurred. The depletion time was noted and the battery was recharged. This procedure was repeated three times.

3.2.6 Electrical safety test

Electrical Safety Test Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 00003850
Military item number: None

Options installed: None

Date of test: 27 Oct 88

Measurements:

Grounding conductor resistance (milliohms): 53.3

Leakage current - Case to ground (microamperes):

unit off, grounded, normal polarity	0.1
unit off, ungrounded, normal polarity	11.2
unit off, ungrounded, reverse polarity	12.1
unit on, grounded, normal polarity	0.1
unit on, ungrounded, normal polarity	11.2
unit on, ungrounded, reverse polarity	12.1

MAXIMUM LIMITS:

ground resistance (milliohms):	150
current (microamperes)	
current (grounded, type A unit):	10
current (ungrounded, type A unit):	100
current (grounded, type B unit):	50
current (ungrounded, type B unit):	500

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.7 Human factors evaluation

Human Factors Evaluation Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 0000385
Military item number: None

Options installed: None

Date of test: 27 Oct 88

Item configuration during test: Item prepared for operation,
sitting on a countertop.

Checklist for HFE

RESULTS

VISUAL DISPLAYS:

Satisfactory

display type, format, content
location of displays
indicator lights
scalar displays
color coding
legends and labels
cathode ray tubes
counters
flags, go/no go, center-null indicators

Comments: None

CONTROLS:

Satisfactory

location
characteristics of controls
labeling
control - display relationships

Comments: None

TIME REQUIRED TO PREPARE FOR OPERATION (list in comment)

Comments: Less than 5 minutes.

MAINTAINABILITY:

Satisfactory

- component location
- component characteristics
- rests and stands
- covers, cases, access doors
- handles
- lubrication
- component mounting
- cord storage provisions
- external accessibility
- internal accessibility
- list special tools required
- list realistic inspection requirements
- list realistic inspection intervals

Comments: No external calibration adjustments

CONDUCTORS:

Satisfactory

- binding and securing
- length
- protection
- routing
- conductor coding
- fabrication
- connectors

Comments: None

FASTENERS:

Satisfactory

- access through inspection panel covers
- enclosure fasteners
- device mounting bolts and fasteners

Comments: None

TEST POINTS:

n/a

- general
- location and mounting
- test point labeling and coding

Comments: None

TEST EQUIPMENT:

Satisfactory

general
equipment self-test
indicators (list in comments)
controls
positive indication of proper operation

Comments: Test "T" connector not provided with unit

FUSES AND CIRCUIT BREAKERS:

Unsatisfactory

external accessibility
easy replacement or reset by operator

Comments: No fuses are externally accessible.

LABELS AND CODING:

Satisfactory

placed above controls and displays
near or on the items they identify
not obscured by other equipment components
describe the function of the items they identify
readable from normal operating distance
conspicuous placards adjacent to hazardous items

Comments: None

SAFETY:

Satisfactory

manual
materials
fire and explosive protection
operator protection from mechanical hazards
patient protection from mechanical hazards
electrical safety (operator and patient)

Comments: None

3.2.8 Altitude test

Altitude Test Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 00003850
Military item number: None

Options installed: None

Date of test: 25 Oct 88

Item configuration during test: Item turned on in the standby mode, operating on DC (battery) power, sitting on chamber floor.

Performance test criteria: Consistent and accurate displays and measurements

Ambient conditions outside chamber:

Temperature	71 degrees F
Humidity	57 percent RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None (battery)
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	Serial Port

IN-TEST DATA

Time of test start: 0908

POSTTEST DATA

Posttest performance check (complete check of item and accessories):

Time of test end : 1035

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): tested with
test and evaluation item 7.

Comments on other data: None

3.2.9 Vibration test

Vibration Test Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 00003850
Military item number: None

Options installed: None

Date of test: 24 Oct 88

Item configuration during test: Item strapped down on vibration table fixture; AC and DC operation.

Performance test criteria: Consistent and accurate measurements and displays.

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 VAC
list connections to simulators	None
list connections to dummy loads	Test engineer's arm
list unconnected terminals	Serial port

Ambient conditions

Temperature	73 degrees F
Humidity	76 percent RH
Barometric pressure	1 atm

IN-TEST DATA

Data and performance checks during test:

Times and dates of test start: 0800

Time at first check:

X: 0800 Y: 0905 Z: 1325

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Time at second check:

X: 0855

Y: 1000

Z: 1430

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

POSTTEST DATA

Posttest performance check (complete check of item and accessories):

Item functional (based on performance test criteria): Yes

Item intact: Yes

Deviation from pretest: None

Comments on item setup or checks:

Times are on different days

Comments on test run (including interruptions):

Z-axis test aborted at 14 minutes because a strap loosened. Test was restarted after strap was tightened.

Comments on other data: None

3.2.10 High temperature test

High Temperature Test (Equipment Operating) Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 00003850
Military item number: None

Options installed: None

Date of test: 25 Nov 88

Item configuration during test: Unit was sitting on chamber floor, ready for operation.

Performance test criteria: Consistent and accurate displays and measurements.

Ambient conditions outside chamber:

Temperature	23 degrees C
Humidity	63 percent RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check :

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 VAC
list connections to simulators	None
list connections to dummy loads	Test engineer's arm
list unconnected terminals	Serial port
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.6
distance from floor (meters)	0.0

IN-TEST DATA

Time of test start: 0815

Performance checks during test:

First check:

Time: 06
Temperature: 49 degrees C
Humidity: 16 percent RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria):
Yes, all OK
Deviation from pretest: None

Second check:

Time: 0930
Temperature: 49 degrees C
Humidity: 16 percent RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria):
Yes, all OK
Deviation from pretest: None

Third check:

Time: 1000
Temperature: 49 degrees C
Humidity: 15 percent RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria):
Yes, all OK
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1100
Item functional (based on performance test criteria):
Yes, all OK
Deviation from pretest: None

Comments on item setup or checks:

tested at same time with test and evaluation program
item 7.

Comments on test run (including interruptions): None

Comments on other data: None

3.2.11 High temperature storage test

High Temperature Test (Equipment in Storage) Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 00003850
Military item number: None

Options installed: None

Date of test: 29 Nov 88

Item configuration during test: Sitting on chamber floor, in storage, not operating.

Performance test criteria: Consistent and accurate displays and measurements.

Ambient conditions outside chamber:

Temperature	20 degrees C
Humidity	43 percent RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	all
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.6
distance from floor (meters)	0.0

Time of test start: 0820

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)

Time of test end: 1445
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data:
Tested at same time with test and evaluation program
item 7.

3.2.12 Low temperature test

Low Temperature Test (Equipment Operating) Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 00003850
Military item number: None

Options installed: None

Date of test: 25 Nov 88

Item configuration during test: Sitting on chamber floor, ready for operation.

Performance test criteria: Consistent and accurate displays and measurements.

Ambient conditions outside chamber:

Temperature	23 degrees C
Humidity	20 percent RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes
All OK Pass

Installation of item in test facility:

list connections to power	120 VAC
list connections to simulators	None
list connections to dummy loads	Test engineer's arm
list unconnected terminals	Serial port
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.0
distance from floor (meters)	0.0

Time of test start: 1200

Performance checks during test:

First check:

Time: 1230
Temperature: 0 degrees C
Humidity: n/a
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Second check:

Time: 1300
Temperature: 0 degrees C
Humidity: n/a
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Third check:

Time: 1330
Temperature: 0 degrees C
Humidity: n/a
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)

Time of test end: 1405
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Comments on item setup or checks:

Tested at same time with test and evaluation program item 7.

Comments on test run (including interruptions):

Condensation on EUT was allowed to dry before final performance check.

Comments on other data: None

3.2.13 Low temperature storage test

Low Temperature Test (Equipment in Storage) Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 0000385
Military item number: None

Options installed: None

Date of test: 1 Dec 88

Item configuration during test: AC power cord and cuff tube coiled and laying on top of the unit. The unit is in storage, not operating.

Performance test criteria: Consistent and accurate displays and measurements

Ambient conditions outside chamber:

Temperature	20 degrees C
Humidity	49 percent RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	All
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.6
distance from floor (meters)	0.0

Time of test start: 0807

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1434

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: The unit was allowed to return to ambient conditions overnight before final performance check.

3.2.14 Humidity test

Humidity Test Report Form

Nomenclature: Blood pressure monitor
Manufacturer: Physio Control
Model number: LIFESTAT® 200
Serial number: 00003850
Military item number: None

Options installed: None

Date of test: 5 Dec 88

Item configuration during test: The unit was sitting on the chamber floor, ready for operation.

Performance test criteria: Consistent and accurate displays and measurements.

Ambient conditions outside chamber:

Temperature	19 degrees C
Humidity	49 percent RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 VAC
list connections to simulators	Test engineer's arm
list connections to dummy loads	None
list unconnected terminals	None
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.6
distance from floor (meters)	0.0

IN-TEST DATA

Time of test start: 1110

Performance checks during test:

First check:

Time: 1155
Temperature: 29.5 degrees C
Humidity: 95 percent RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Second check:

Time: 1240
Temperature: 29.5 degrees C
Humidity: 95 percent RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Third check:

Time: 1325
Temperature: 29.5 degrees C
Humidity: 95 percent RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Fourth check:

Time: 1410
Temperature: 29.5 degrees C
Humidity: 95 percent RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Fifth check:

Time: 1455
Temperature: 29.5 degrees C
Humidity: 95 percent RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1550

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data:

This test item was tested with test and evaluation item 7.

3.2.15 Electromagnetic characteristics test

Electromagnetic Characteristics Testing Evaluation of Performance

T & E Item Number: 06

Date: 1 Nov 1988

Nomenclature: Blood pressure monitor

Manufacturer: Physio Control

Model number: LIFESTAT® 200

Serial number: 00003850

Military item number: n/a

Conducted emissions tests

CE01 Testing configuration(s): n/a
 Performance (pass/fail): n/a

 Comments: n/a

CE02 Testing configuration(s): Operating on copper work
 bench.
 Performance (pass/fail): Pass

 Comments: Both hot and neutral conductors tested.

CE04 Testing configuration(s): Operating on copper work
 bench.
 Performance (pass/fail): Fail

 Comments: NB failure 0.2 to 13.3 dB over specifi-
 cations in range 13 to 53 MHz; BB failure 3 to
 4.5 dB over specifications in range 15.6 to 16 MHz.

Conducted susceptibility tests

CS02 Testing configuration(s): Operating on test bench,
 connected to test jig.
 Performance (pass/fail): n/a

 Comments: Unable to test because noise generated
 by the unit is greater than the test signal (unable
 to measure test signal).

Radiated emissions tests

Radiated susceptibility tests

3-35

3.3 CRITERIA, SIGNIFICANT PROBLEMS, AND SUGGESTED IMPROVEMENTS

3.3.1 Criteria

Item			<u>Applicable</u>
<u>No.</u>	<u>Criteria (Source)</u>	<u>Remarks</u>	<u>subparagraph</u>
1	The physical inventory is conducted solely for investigation and documentation.	N/A	2.1.2.1
2	The LIFESTAT® 200 will display consistent and accurate measurements.	met	2.1.2.2
3	Verify manufacturer's specified full power internal battery life expectancy of 30 minutes during continuous operation in the 2-minute cycle.	met	2.2.2
4	The LIFESTAT® 200 will meet the limits established in NAEP 99 for electrical safety of medical equipment.	met	2.3.2
5	The LIFESTAT® 200 will be rated satisfactory in all major categories of the evaluation. These include: Visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.	partially met	2.4.2
6	The LIFESTAT® 200 will display consistent and accurate measurements while exposed to an altitude equivalency of 15,000 feet above sea level.	met	2.5.2
7	The LIFESTAT® 200 will remain operational and display consistent and accurate measurements while exposed to vibrational stresses.	met	2.6.2

8	The LIFESTAT® 200 will display consistent and accurate measurements during the high temperature operation check.	met	2.7.2.1
9	The LIFESTAT® 200 will display consistent and accurate measurements after the high temperature storage.	met	2.7.2.2
10	The LIFESTAT® 200 will display consistent and accurate measurements during the low temperature operation check.	met	2.8.2.1
11	The LIFESTAT® 200 will display consistent and accurate measurements after the low temperature storage.	met	2.8.2.2
12	The LIFESTAT® 200 will display consistent and accurate measurements while exposed to a high humidity.	met	2.9.2
13	The LIFESTAT® 200 will not produce emissions in excess of the limits set forth in MIL-STD-461A Notice 4, paragraph 6.13.	partially met	2.10.2.1

14	The LIFESTAT® 200 will not malfunction when it is subjected to radiated fields as specified in MIL-STD-461A, Notice 4, paragraph 6.20.	met	2.10.2.2
15	The LIFESTAT® 200 will not conduct emissions in excess of the limits set forth in paragraphs 6.1 and 6.2, MIL-STD-461A, Notice 4.	partially met	2.10.2.3
16	The LIFESTAT® 200 will not malfunction when it is subjected to conducted emissions as specified in MIL-STD-461A, Notice 4, paragraphs 6.7 and 6.10.	partially met	2.10.2.4
17	The flight surgeon will be able to operate the LIFESTAT® 200 without physical or functional restrictions aboard the aircraft.	partially met	2.11.2.1
18	The LIFESTAT® 200 will not radiate EMI to disrupt or interfere with the other equipment or systems aboard the aircraft.	met	2.12.2.2
19	The aircraft will not radiate EMI to disrupt or interfere with the LIFESTAT® 200.	met	2.12.2.3

3.3.2 Significant problems which require corrective action

None

3.3.3 Suggested improvements

None

3.4 REFERENCES

- 3.4.1 Department of Defense. 1971. EMI characteristics, requirements for equipment. Washington, D.C. MIL-STD-461A, Notice 4. February.
- 3.4.2 Department of Defense. 1971. EMI characteristics, measurement of. Washington, D.C. MIL-STD-462, Notice 3. February.
- 3.4.3 Department of Defense. 1983. Environmental test methods and engineering guidelines. Washington, D.C. MIL-STD-810D. July.
- 3.4.4 Department of the Army. 1982. Environmental protection and enhancement. Washington, D.C. AR-200-1. June.
- 3.4.5 Department of the Army. 1987. Maintenance management procedures for medical equipment. Washington, D.C. TB 38-750-2. April.
- 3.4.6 Department of Defense. 1985. Standard general requirements for electronic equipment. Washington, D.C. MIL-STD-454K. February.
- 3.4.7 Underwriters Laboratory's, Inc. 1978. Standard for safety, medical and dental equipment. Chicago, Illinois. UL-544.
- 3.4.8 Department of Defense. 1989. Human engineering design criteria for military systems, equipment, and facilities. Washington, D.C. MIL-STD-1472D. March.
- 3.4.9 Association for the Advancement of Medical Instruments. Human factors engineering guidelines and preferred practices for the design of medical devices. Arlington, Virginia. AAMI-HE-1988. February.
- 3.4.10 Department of the Army. 1978. Operator's manual, UH-60 and EH-60 helicopter, with changes 1-5. Washington, D.C. TM 55-1520-237-10. January.
- 3.4.11 Department of the Army. 1987. Maintenance management procedures for medical equipment. Washington, D.C. TB 38-750-2. April.
- 3.4.12 National Fire Protection Association. 1987. Standard for health care facilities. Quincy, Massachusetts. February.

3.4.13 Physio Control. 1987. Operating instructions, LIFESTAT® 200 noninvasive blood pressure monitor. Redmond, Washington. P/N 802611-02.

3.4.14 Mitchell, G. W., and Adams, J. E. 1988. Technical test and evaluation of aeromedical equipment. Fort Rucker, AL: U.S. Army Aeromedical Research Laboratory. USAARL Letter Report LR-88-16-1-2.

3.5 ABBREVIATIONS

AVSCOM	Army Aviation Systems Command
AEST	aeromedical equipment suitability test
AWR	airworthiness release
BB	broadband
BPM	beats per minute
CAAF	Cairns Army Airfield
CRT	cathode ray tube
dB	decibel
DC	direct current
ECG	electrocardiograph
EMC	electromagnetic compatibility
EMI	electromagnetic interference
EUT	equipment under test
fpm	feet per minute
GFE	government furnished equipment
Gpk	gravity, peak
G(rms)	gravity (root mean square)
Hz	hertz
IAW	in accordance with
ITOP	in-flight test operating procedure
IGE	in-ground effect
kHz	kilohertz
KIAS	knots indicated airspeed
LCD	liquid crystal display
LED	light emitting diode
LIFESTAT® 200	Physio Control blood pressure monitor, model LIFESTAT® 200
MEDEVAC	medical evacuation
MHz	mega hertz
MIL-STD	military standard
ml	milliliter
mm	millimeter
mmHg	millimeters of Mercury
MSL	mean sea level
NAFP	National Association of Fire Prevention
NB	narrowband
NBC	nuclear, biological and chemical

NVG	night vision goggle
RAM	random access memory
RF	radio frequency
RH	relative humidity
ROM	read only memory
TB	technical bulletin
TFT	technical feasibility testing
T & E	test and evaluation
UES	Universal Energy Systems, Inc.
USAARL	U.S. Army Aeromedical Research Laboratory
V/m	volts per meter

3.6 LIST OF MANUFACTURERS

- 3.6.1 Physio-Control Corporation
 11811 Willows Road Northeast
 Post Office Box 97006
 Redmond, WA 98073-9706
- 3.6.2 Sikorsky Aircraft
 6900 Main Street
 Stratford, CT 06601
- 3.6.3 Neurodyne-Dempsey,, Inc
 200 Arrowhead Drive
 Carson City, NV 89701
- 3.6.4 Tenney Engineering, Inc.
 1090 Springfield Road
 Post Office Box 3142
 Union, NJ 07083
- 3.6.5 Uholtz-Dickey Corporation
 6 Brookside Drive
 Wallingford, CT 06492
- 3.6.6 BioTek Instruments, Inc.
 1 Mill Street
 Burlington, VT 05401

3.7 DISTRIBUTION LIST

Commander, U.S. Army Natick Research,
Development and Evaluation Center
ATTN: STRNC-MIL (Documents
Librarian)
Natick, MA 01760-5040

Naval Submarine Medical
Research Laboratory
Medical Library, Naval Sub Base
Box 900
Groton, CT 06340

Commander/Director
U.S. Army Combat Surveillance
and Target Acquisition Lab
ATTN: DELCS-D
Fort Monmouth, NJ 07703-5304

Commander
10th Medical Laboratory
ATTN: Audiologist
APO New York 09180

Naval Air Development Center
Technical Information Division
Technical Support Detachment
Warminster, PA 18974

Commanding Officer, Naval Medical
Research and Development Command
National Naval Medical Center
Bethesda, MD 20814-5044

Deputy Director, Defense Research
and Engineering
ATTN: Military Assistant
for Medical and Life Sciences
Washington, DC 20301-3080

Commander, U.S. Army Research
Institute of Environmental Medicine
Natick, MA 01760

U.S. Army Avionics Research
and Development Activity
ATTN: SAVAA-P-TP
Fort Monmouth, NJ 07703-5401

U.S. Army Communications-Electronics
Command
ATTN: AMSEL-RD-ESA-D
Fort Monmouth, NJ 07703

Library
Naval Submarine Medical Research Lab
Box 900, Naval Sub Base
Groton, CT 06349-5900

Commander
Man-Machine Integration System
Code 602
Naval Air Development Center
Warminster, PA 18974

Commander
Naval Air Development Center
ATTN: Code 602-B (Mr. Brindle)
Warminster, PA 18974

Commanding Officer
Harry G. Armstrong Aerospace
Medical Research Laboratory
Wright-Patterson
Air Force Base, OH 45433

Director
Army Audiology and Speech Center
Walter Reed Army Medical Center
Washington, DC 20307-5001

Commander, U.S. Army Institute
of Dental Research
ATTN: Jean A. Setterstrom, Ph. D.
Walter Reed Army Medical Center
Washington, DC 20307-5300

Naval Air Systems Command
Technical Air Library 950D
Room 278, Jefferson Plaza II
Department of the Navy
Washington, DC 20361

Naval Research Laboratory Library
Shock and Vibration
Information Center, Code 5804
Washington, DC 20375

Director, U.S. Army Human
Engineering Laboratory
ATTN: Technical Library
Aberdeen Proving Ground, MD 21005

Commander, U.S. Army Test
and Evaluation Command
ATTN: AMSTE-AD-H
Aberdeen Proving Ground, MD 21005

Director
U.S. Army Ballistic
Research Laboratory
ATTN: DRXBR-OD-ST Tech Reports
Aberdeen Proving Ground, MD 21005

Commander
U.S. Army Medical Research
Institute of Chemical Defense
ATTN: SGRD-UV-AO
Aberdeen Proving Ground,
MD 21010-5425

Commander, U.S. Army Medical
Research and Development Command
ATTN: SGRD-RMS (Ms. Madigan)
Fort Detrick, Frederick, MD 21702-5012

Director
Walter Reed Army Institute of Research
Washington, DC 20307-5100

HQ DA (DASG-PSP-O)
5109 Leesburg Pike
Falls Church, VA 22041-3258

Naval Research Laboratory
Library Code 1433
Washington, DC 20375

Harry Diamond Laboratories
ATTN: Technical Information Branch
2800 Powder Mill Road
Adelphi, MD 20783-1197

U.S. Army Materiel Systems
Analysis Agency
ATTN: AMXSY-PA (Reports Processing)
Aberdeen Proving Ground
MD 21005-5071

U.S. Army Ordnance Center
and School Library
Simpson Hall, Building 3071
Aberdeen Proving Ground, MD 21005

U.S. Army Environmental
Hygiene Agency
Building E2100
Aberdeen Proving Ground, MD 21010

Technical Library Chemical Research
and Development Center
Aberdeen Proving Ground, MD
21010-5423

Commander
U.S. Army Medical Research
Institute of Infectious Disease
SGRD-UIZ-C
Fort Detrick, Frederick, MD 21702

Director, Biological
Sciences Division
Office of Naval Research
600 North Quincy Street
Arlington, VA 22217

Commander
U.S. Army Materiel Command
ATTN: AMCDE-XS
5001 Eisenhower Avenue
Alexandria, VA 22333

Commandant
U.S. Army Aviation
Logistics School ATTN: ATSQ-TDN
Fort Eustis, VA 23604

Headquarters (ATMD)
U.S. Army Training
and Doctrine Command
Fort Monroe, VA 23651

Structures Laboratory Library
USARTL-AVSCOM
NASA Langley Research Center
Mail Stop 266
Hampton, VA 23665

Naval Aerospace Medical
Institute Library
Building 1953, Code 03L
Pensacola, FL 32508-5600

Command Surgeon
HQ USCENTCOM (CCSG)
U.S. Central Command
MacDill Air Force Base FL 33608

Air University Library
(AUL/LSE)
Maxwell Air Force Base, AL 36112

U.S. Air Force Institute
of Technology (AFIT/LDEE)
Building 640, Area B
Wright-Patterson
Air Force Base, OH 45433

Henry L. Taylor
Director, Institute of Aviation
University of Illinois-Willard Airport
Savoy, IL 61874

COL Craig L. Urbauer, Chief
Office of Army Surgeon General
National Guard Bureau
Washington, DC 50310-2500

Commander
U.S. Army Aviation Systems Command
ATTN: SGRD-UAX-AL (MAJ Gillette)
4300 Goodfellow Blvd., Building 105
St. Louis, MO 63120

U.S. Army Aviation Systems Command
Library and Information Center Branch
ATTN: AMSAV-DIL
4300 Goodfellow Boulevard
St. Louis, MO 63120

Federal Aviation Administration
Civil Aeromedical Institute
Library AAM-400A
P.O. Box 25082
Oklahoma City, OK 73125

Commander
U.S. Army Academy
of Health Sciences
ATTN: Library
Fort Sam Houston, TX 78234

Commander
U.S. Army Institute of Surgical Research
ATTN: SGRD-USM (Jan Duke)
Fort Sam Houston, TX 78234-6200

AAMRL/HEX
Wright-Patterson
Air Force Base, OH 45433

John A. Dellinger,
Southwest Research Institute
P. O. Box 28510
San Antonio, TX 78284

Product Manager
Aviation Life Support Equipment
ATTN: AMCPM-ALSE
4300 Goodfellow Boulevard
St. Louis, MO 63120-1798

Commander
U.S. Army Aviation
Systems Command
ATTN: AMSAV-ED
4300 Goodfellow Boulevard
St. Louis, MO 63120

Commanding Officer
Naval Biodynamics Laboratory
P.O. Box 24907
New Orleans, LA 70189-0407

Assistant Commandant
U.S. Army Field Artillery School
ATTN: Morris Swott Technical Library
Fort Sill, OK 73503-0312

Commander
U.S. Army Health Services Command
ATTN: HSOP-SO
Fort Sam Houston, TX 78234-6000

Director of Professional Services
HQ USAF/SGDT
Bolling Air Force Base, DC 20332-6188

U.S. Army Dugway Proving Ground
Technical Library, Building 5330
Dugway, UT 84022

U.S. Army Yuma Proving Ground
Technical Library
Yuma, AZ 85364

AFFTC Technical Library
6510 TW/TSTL
Edwards Air Force Base,
CA 93523-5000

Commander
Code 3431
Naval Weapons Center
China Lake, CA 93555

Aeromechanics Laboratory
U.S. Army Research and Technical Labs
Ames Research Center, M/S 215-1
Moffett Field, CA 94035

Sixth U.S. Army
ATTN: SMA
Presidio of San Francisco, CA 94129

Commander
U.S. Army Aeromedical Center
Fort Rucker, AL 36362

U.S. Air Force School
of Aerospace Medicine
Strughold Aeromedical Library Technical
Reports Section (TSKD)
Brooks Air Force Base, TX 78235-5301

Dr. Diane Damos
Department of Human Factors
ISSM, USC
Los Angeles, CA 90089-0021

U.S. Army White Sands
Missile Range
ATTN: STEWS-IM-ST
White Sands Missile Range, NM 88002

U.S. Army Aviation Engineering
Flight Activity
ATTN: SAVTE-M (Tech Lib) Stop 217
Edwards Air Force Base, CA 93523-5000

Ms. Sandra G. Hart
Ames Research Center
MS 262-3
Moffett Field, CA 94035

Commander, Letterman Army Institute
of Research
ATTN: Medical Research Library
Presidio of San Francisco, CA 94129

COL Eugene S. Channing, O.D.
Brooke Army Medical Center
ATTN: HSHE-EAH-O
Fort Sam Houston, TX 78234-6200

Commander
U.S. Army Medical Materiel
Development Activity
Fort Detrick, Frederick, MD 21702-5009

Commander
U.S. Army Aviation Center
Directorate of Combat Developments
Building 507
Fort Rucker, AL 36362

U. S. Army Research Institute
Aviation R&D Activity
ATTN: PERI-IR
Fort Rucker, AL 36362

Commander
U.S. Army Safety Center
Fort Rucker, AL 36362

U.S. Army Aircraft Development
Test Activity
ATTN: STEBG-MP-P
Cairns Army Air Field
Fort Rucker, AL 36362

Commander U.S. Army Medical Research
and Development Command
ATTN: SGRD-PLC (COL Sedge)
Fort Detrick, Frederick, MD 21702

MAJ John Wilson
TRADOC Aviation LO
Embassy of the United States
APO New York 09777

Netherlands Army Liaison Office
Building 602
Fort Rucker, AL 36362

British Army Liaison Office
Building 602
Fort Rucker, AL 36362

Italian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Directorate of Training Development
Building 502
Fort Rucker, AL 36362

Chief
USAHEL/USAAVNC Field Office
P. O. Box 716
Fort Rucker, AL 36362-5349

Commander U.S. Army Aviation Center
and Fort Rucker
ATTN: ATZQ-CG
Fort Rucker, AL 36362

Commander/President
TEXCOM Aviation Board
Cairns Army Air Field
Fort Rucker, AL 36362

MAJ Terry Newman
Canadian Army Liaison Office
Building 602
Fort Rucker, AL 36362

German Army Liaison Office
Building 602
Fort Rucker, AL 36362

LTC Patrice Cottebrune
French Army Liaison Office
USAAVNC (Building 602)
Fort Rucker, AL 36362-5021

Brazilian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Australian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Dr. Garrison Rapmund
6 Burning Tree Court
Bethesda, MD 20817

Commandant Royal Air Force
Institute of Aviation Medicine
Farnborough Hants UK GU14 6SZ
Dr. A. Kornfield, President
Biosearch Company
3016 Revere Road
Drexel Hill, PA 29026

Commander
U.S. Army Biomedical Research
and Development Laboratory
ATTN: SGRD-UBZ-I
Fort Detrick, Frederick, MD 21702

Defense Technical Information Center
Cameron Station
Alexandra, VA 22313

Commander, U.S. Army Foreign Science
and Technology Center
AIFRTA (Davis)
220 7th Street, NE
Charlottesville, VA 22901-5396

Director,
Applied Technology Laboratory
USARTL-AVSCOM
ATTN: Library, Building 401
Fort Eustis, VA 23604

U.S. Army Training
and Doctrine Command
ATTN: Surgeon
Fort Monroe, VA 23651-5000

Aviation Medicine Clinic
TMC #22, SAAF
Fort Bragg, NC 28305

U.S. Air Force Armament
Development and Test Center
Eglin Air Force Base, FL 32542

Commander, U.S. Army Missile
Command
Redstone Scientific Information Center
ATTN: AMSMI-RD-CS-R/ILL
Documents Redstone Arsenal, AL 35898

U.S. Army Research and Technology
Laboratories (AVSCOM)
Propulsion Laboratory MS 302-2
NASA Lewis Research Center
Cleveland, OH 44135

Dr. H. Dix Christensen
Bio Medical Science Building, Room 753
Post Office Box 26901
Oklahoma City, OK 73190

Dr. Christine Schlichting
Behavioral Sciences Department
Box 900, NAVUBASE NLON
Groton, CT 06349-5900

Commander
U.S. Army Aviation Systems Command
ATTN: AMSAV-ECU
4300 Goodfellow Bouvelard
St. Louis, MO 63120-1790

Commandant
Academy of Health Sciences
ATTN: HSHA-COM (LTC Huether)
Fort Sam Houston, TX 78234

U.S. Air Force Armament
Development and Test Center
Eglin Air Force Base, FL 32542

COL Eugene S. Channing, O.D.
Brooke Army Medical Center
ATTN: HSHE-EAH-O
Fort Sam Houston, TX 78234-6200